
# NIKKEN®

BIAXIAL  
**Body Energizer™**  
(Énergisant Physique Biaxial<sup>MC</sup>)

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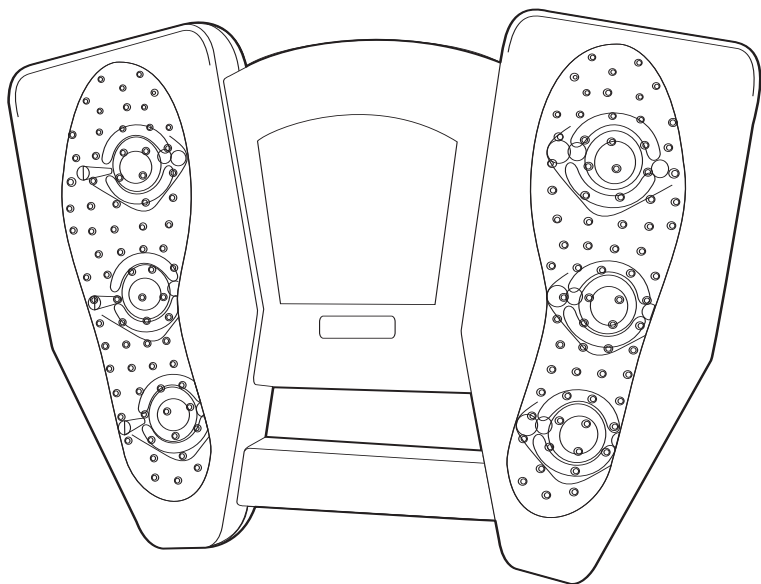## **INSTRUCTION MANUAL**

When you put your feet on the Nikken Biaxial Body Energizer, you can begin to feel a difference almost immediately.

The gentle, soothing vibration and the specially textured surface provide a mild massage effect. Within minutes this soothing sensation begins to spread. Soon, this sense of relaxation and release of tension suffuses through your body.

The Biaxial Body Energizer incorporates a patented magnetic technology available only from Nikken — Magnetic Biaxial Rotation.

Magnetic Biaxial Rotation uses a revolving magnet, that also rotates on a second axis. The magnetic field constantly changes its angle, as much as 25 times per second. This creates three-dimensional movement of the field. No other magnetic device can achieve this effect.

The Nikken Biaxial Body Energizer contains six complete Magnetic Biaxial Rotation units, three under each foot. Each comprises a strong 1,800-gauss magnetic sphere and an individual motor assembly.

The simple and rugged direct-drive motor system minimizes moving parts, to make it quiet and reliable.

The Biaxial Body Energizer features selectable speed from 500 to 1,500 rpm. A top-mounted LED readout shows the speed setting.


The foot surfaces are textured to provide a feeling of comfort and gentle, stimulating massage. These surfaces are slightly inclined at an ergonomically precise angle.

The durable case is made of high-impact Lexan®. Built into the case is a convenient carrying handle.

The Biaxial Body Energizer operates using household AC current. A plug-in adapter is included.

The adapter is similar to that provided with the Nikken PalmMag™ 1000. Either adapter may be used with both products.

## OPERATING THE BIAXIAL BODY ENERGIZER

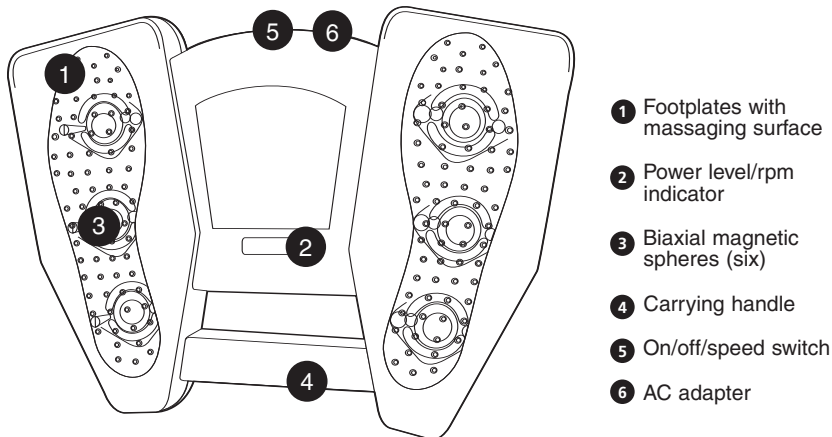You can use the Biaxial Body Energizer while working at a desk, watching television, relaxing — any activity when you want to put up your feet and be comfortable.

The Biaxial Body Energizer is designed for use while seated. **DO NOT** stand on the unit when it is operating or when it is switched off. Do not sit on it.

**DO NOT** place excessive downward force on the unit while it is operating. Remove feet from the Biaxial Body Energizer when turning it on or off or when adjusting the speed setting. Leaning forward to reach these controls with feet on the unit may cause the top surface to come into contact with the rotating spheres and stop their motion.

Place the Biaxial Body Energizer in front of a chair or seat. Make sure the energizer is placed so that the front of the unit points forward.

Connect the adapter plug to a standard 100/120-volt outlet. Insert the other end of the adapter cord into the receptacle on the unit.

For best results, remove your shoes before placing your feet on the Body Energizer. Socks or stockings may be worn if desired.

Locate the power switch and turn it to activate the unit. The speed is selectable from 500 to 1,500 revolutions per minute. An LED readout on the top surface indicates the speed when operating.

After some experience with the Biaxial Body Energizer, you can find the power/speed that is most comfortable for you.

Nikken recommends using the Energizer for a period of 10-15 minutes with the unit activated. You can also keep your feet on the Biaxial Body Energizer with the power turned off.

## CARE INSTRUCTIONS

Disconnect the Biaxial Body Energizer from the power source when not in use. Store the Energizer and adapter in a cool, dry place.

To clean the Energizer, wipe gently with a damp cloth. The unit may also be wiped with a mild disinfectant solution. Do not immerse the Energizer in water or subject it to a stream of water. Do not use chemical solvents or any materials that can damage the case.

Do not attempt to open the Energizer case. There are no user serviceable components inside.

Do not place weight on the Energizer. Use only while seated. Do not lean forward when the Biaxial Body Energizer is in operation.

## SPECIFICATIONS

### Magnetic units

6, each 1,800 gauss, 6-volt motors

### Speed

Selectable, 500-1,500 rpm

### Weight

3 lb/1.36 kg

## TROUBLESHOOTING

If difficulty is encountered when attempting to use the Biaxial Body Energizer, here are areas for attention.

**Situation:** Unit will not operate when switch is turned on.

**Remedy:** Make sure that adapter is firmly connected at both ends and check that the electric outlet is powered.

**Situation:** Unit is operating but one or more spheres are not rotating, or not rotating in both directions.

**Remedy:** Turn off the Biaxial Body Energizer for a few seconds to reset it, then restart at slow speed and increase speed gradually.

## WARNING

Do not use the Biaxial Body Energizer if you have a pacemaker or other sensitive implanted device. Women in the first trimester of pregnancy should consult a physician before using the Biaxial Body Energizer. Use of the Biaxial Body Energizer may affect other medical instruments or devices. Placing the Biaxial Body Energizer in direct contact with some watches, audio/video tapes, credit cards or electronic equipment may impair their use.

Operate the Biaxial Body Energizer only using the approved Body Energizer or PalmMag 1000 AC adapter. Use of other adapters or devices may cause damage to the Biaxial Body Energizer or create a fire hazard.

## WARRANTY

The Biaxial Body Energizer is warranted against manufacturing defects for 90 days after date of purchase.

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©2005 Nikken Inc. • Printed in USA / Imprimé aux É.-U.
LEXAN is a registered trademark of GE Plastics, General Electric Corporation.
Nikken Inc.
Nihon Kenko Zoushin Kenkyukai, Canada Corp.